

Instruction concerning orthopaedic Implants (Bone Plates & Bone Screws) made by Miraclus Orthotech Pvt Ltd., located at Survey No 106, Village Vavdi, Opp Safex Industries, Behind Elite Boarding School, Vavdi, Ta. Kheda-387550, Gujrat, India

CONTENTS

The Device package contains single use implant (Bone Plate& Bone Screws) of the Miraclus Orthotech Pvt. Ltd.

DESCRIPTION

The Bone plates & Bone Screws are single use device supplied Non-sterile. The devices are available in different sizes.

INTENDED USE

Bone Plates:

The bone plates is intended for internal fixation of fractures and reconstruction of bones including the Scapula, Olecranon, Humerus, Radius, Ulna, Pelvis, Distal, Tibia, Fibula, Hand and Foot in adults and for long bones in adolescents in whom the growth plates are fused. Examples of these internal fixations and reconstructions include compression fractures, intra-articular and extra-articular fractures, displaced fractures, osteotomies, non-unions and mal-unions. This system can be used for ventral, dorsal or orthogonal application.

Bone Screws:

Bone Screws is intended for internal fixation of fractures and reconstruction of bones including the Humerus, Tibia, Fibula, Hand and Foot in adults and for long bones in adolescents. Bone Screws fasten plates to bone, maintain bone fragments in their relative position, or hold together fragments of bone (in the case of lag screws). They vary according to the manner in which they are inserted into bone, their function, size, the type of bone they are intended for, and the manner in which they couple with the screw driver (i.e. slot, cruciate, square, hexagon).

Correct selection of the implants is extremely important:

- Responsibility of the proper selection of patients, adequate training, experience in the choice, placement of the implant & the decision to leave or remove implant postoperatively, rests with the surgeon.
- Our Bone Plates & Bone Screws are available in variety of configurations, these shall be used in combination with related corresponding implants & instruments made by Miraclus Orthotech Pvt. Ltd. only.
- The product should be used in combination with the devices made up similar material only. (Titanium Gr.23 implants with Titanium Gr.23 or Titanium Grade 5 & SS 316LVM with SS 316LVM).
- The Surgeon should discuss the expectation of the surgery inherent the use of the product with the patient. Particular attention should be given to a discussion postoperatively & the necessity should be focused for periodic medical follow-up.
- The Correct selection of the product is extremely important. The product should be used in the correct anatomical location, consistent with the accepted standard for the internal fixation. Failure to use the appropriate product for the application may result in a premature clinical failure. Failure to use the proper component to ensure adequate blood supply & provide rigid fixation may result in loosening, bending or cracking of the product and/ or bone fracture.

Note: Surgical techniques & combination chart will be available upon doctor request.

CONTRAINDICATIONS

Do not use the Orthopedics implants in cases of:

- Inadequate bone quantity and/or bone quality
- Hypersensitivity to metal or allergic reaction
- Early or Late Inspection, both deep and / or superficial
- Patients with limited blood supply
- Patient within whom co-operation or mental competence is lacking, thereby reducing patient compliance.

ADVERSE REACTIONS

- Adverse reactions may include but are not limited to:
 - Clinical failure (i.e. pain or injury) due to bending, loosening, breakage of implant, loos fixation, dislocation and/or migration
 - Pain, discomfort, and/or abnormal sensations due to the presence of the implant.

- Primary and/or secondary infections.
- Allergic reactions to implant material.
- Necrosis of bone or decrease of bone density.
- Injury to vessels, nerves and organs.
- Elevated fibrotic tissue reaction around the surgical area

SAFETY PRECAUTIONS

- The Product should only be used by the medical personnel who hold relevant qualification.
- Never use the product that has been damaged by Improper handling in the hospital or in any other way.
- Never reuse an implant. Although the implant appears to be undamaged, previous stresses may have created non-visible damage that could result in implant failure.
- Safety Precaution for Special Cases

Pregnant Women

- ✓ Ensure that there should be less blood loss during the surgery.
- ✓ Anaesthesia should not be used in such case.
- ✓ Operational environment must be free from radiation.

Infant / Children

- ✓ Ensure that there should be less blood loss during the surgery.
- ✓ Operational environment must be free from radiation.
- ✓ Epiphysis should not be damaged

HOW SUPPLIED/STORAGE:

- The implants are individually packed in protective packaging that is labelled to its contents properly. All Single use **Non-sterile** implants are supplied.
 - Implants should be stored in the original protective packaging.
 - Store the implants in a dry and dust-free place (standard hospital environment).

INSPECTION:

Before use, inspect the box carefully. Do not use when

- Implants has scratches & damage
- Improper threads with damages
- Prior to surgery check suitability of fixation of this implant with its corresponding implant, and also ensure strength of whole assembly.

OPERATING INSTRUCTIONS

The Miraclus Orthotech Pvt. Ltd. implants should be implanted only with the related corresponding instruments made by Miraclus Orthotech Pvt. Ltd.

Also ensure the availability of same implant as standby.

Surgeon should document the implant details (name, item, number, lot number) in surgery record.

PRE-OPERATIVE

- Keep the instructions for use accessible to all staff.
- The operating surgeon must have a thorough understanding of both, the hands-on and conceptual aspects of the established operating techniques. Proper surgical performance of the implantation is the responsibility of the operating surgeon. The operating surgeon draws up an operation plan specifying and documenting the following:
 - Implant component(s) and their dimensions.
 - Determination of intra-operative orientation points.

The following conditions must be fulfilled prior to application:

- All required implant components are sterilized and readily available.
- All requisite sterile implantation instruments must be available and in working order.
- Highly aseptic operating conditions are present.

Sterilization: All Single use NON STERILE implants and instrument used in the surgery must be cleaned & Sterile prior to use.
Remove plastic packing of implant before cleaning.

Cleaning Procedure:

New products must be carefully cleaned before initial sterilization. Only trained personnel must perform cleaning

Equipment: various sized soft-bristled brushes, lint-free cloths, syringes, pipettes and/or water jet, neutral enzymatic cleaner or neutral detergent with a pH 7.

1. Rinse Bone plates & Bone Screws under running cold tap water for a minimum of two minutes. Use a soft-bristled brush to clean the Bone plates & Bone Screws.
2. Soak Bone plates & Bone Screws in a neutral pH enzymatic cleaner or detergent solution for a minimum of ten minutes. Follow the enzymatic cleaner or detergent manufacturer’s instructions for use for correct exposure time, temperature, water quality, and concentration.
3. Rinse Bone plates & Bone Screws with cold water for a minimum of two minutes. Use a syringe, pipette, or water jet to flush lumens, channels, and other hard to reach areas.
4. Manually clean Bone plates & Bone Screws for a minimum of five minutes in a freshly prepared neutral pH enzymatic cleaner or detergent solution using a soft-bristled brush. Clean Bone plates & Bone Screws under water to prevent aerosolization of contaminants.

Note: Freshly prepared solution is a newly-made, clean solution.

5. Rinse Bone plates & Bone Screws thoroughly with deionized (DI) or purified (PURW) water for a minimum of two minutes. Use a syringe, pipette, or water jet to flush lumens and channels.
6. Visually inspect Bone plates & Bone Screws
7. Perform a final rinse on Bone plate using DI or PURW water.
8. Dry Bone plate using a clean, soft, lint-free cloth or clean compressed air.

Note: Cleaning Agent Information: We used the following cleaning agents during internal processes of these cleaning recommendations. These cleaning agents are not listed in preference to other available cleaning agents which may perform satisfactorily- neutral pH enzymatic detergents (e.g. Prolystica 2X Concentrate Enzymatic Cleaner, Enzol, Endozime, and Neodisher Medizym) and neutral pH detergents (e.g. Prolystica 2X Neutral Detergent).

We are suggesting following parameter for the sterilization;

Method	Temperature	Exposure time	Pressure
Steam (autoclave)	121 Deg C.	15 Minutes*	103421 Pa / 0.1 MPa / 15 psi

Note: Recommended Steam Sterilizer (Autoclave) is Class B.

INTRA-OPERATIVE

- Prior to use, verify the integrity of the implant.
- Modification of the Orthopedic Implant Set is not allowed.
- Bending of the Bone plates & Bone Screws set is possible. When contouring these plates, do not over bend and / or bend back in original shape.
- Ensure sufficient rinsing in-situ for cooling and removing of potential wear material.
- Before locking the screw to the plate, the bone has to be correctly repositioned.

POST-OPERATIVE

- Reiterate preoperative instructions to the patient.
- During the post-operative phase, in addition to mobility it is of vital importance that the physician keeps the patient well informed about post surgical behavioral requirements.
- Ensure that the patient is aware of physical activity restrictions and possible adverse reactions.

WARNING:

The use of implants for surgery other than those for which they are intended may result in damage/ breakage of implants or patient injury.

- The operating surgeon and operating room team must be thoroughly familiar with the operating technique, as well as the range of implants and instruments to be applied. Complete information on these subjects must be readily available at the workplace.
- The operating surgeon must be especially trained in orthopedic surgery, biomechanical principles of the skeleton, and the relevant operating techniques.
- The patient is aware of the risks associated with general surgery, orthopedic surgery, and with general anesthesia.
- The patient has been informed about the advantages and disadvantages of the implant & implantation procedure and about possible alternative treatments.
- The implant can be failed due to excessive load, wear and tear or infection.
- Corrective surgery may be necessary if the implant fails.
- The patient must have his/her physician to carry out follow-up examinations of the implants at regular intervals.
- If device used in joints, kindly inform to patient do not move excessively, it may cause pain or damage surrounding tissue where implant were placed.

REVISION SURGERY / IMPLANT REMOVAL

Metallic implants can be loosen, fracture, migrate, cause pain, or stress shield bone even after a fracture is healed, particularly in young active patients. The surgeon must make the final decision on implant removal if either of these occurs. If there are not any of these complications, we recommend the permanent implantation of this implants because of the risk of re-fracture and the possible complications of an additional operation.

- The surgeon must make the final decision on implant removal if either of these occurs;
 - ✓ Choice of Patient
 - ✓ Doctor's Advice based on the clinical condition of the patient
 - ✓ Deep Wound Infection/Bone Atrophy
 - ✓ Growing Skeleton
 - ✓ Tenosynovitis
 - ✓ Intra-Articular Material
 - ✓ Post – traumatic Arthritis
 - ✓ Avascular Necrosis
 - ✓ Intractable Pain
 - ✓ Perforating Material
 - ✓ Infection
 - ✓ Paresthesia
- Time of removal of implant shall be suggested by the doctor depending upon the clinical condition of the patient either after the surgery or during the follow ups.
- Removal of Implant may cause the risk of re-fracture, neurovascular injury & infection.
- Bone in-growth and wear of the implant can make the removal difficult.

MRI SAFETY INFORMATION

Miraclus Orthotech implants are manufactured from SS 316LVM, Titanium Grade 23 & Titanium (Gr-5) material, hence it do not pose any safety risk.

- Patients should be directed to seek a medical opinion before entering potentially adverse environments that could affect the performance of the implants, such as electromagnetic or magnetic fields, including a magnetic fields, including a magnetic resonance environment.
- The Miraclus Orthotech Pvt. Ltd. implants has not been evaluated for safety and compatibility in the MR environment but on the basis of literature study below mentioned points can be taken care during MRI
 - ✓ The minimum recommended time after the implantation that allows patients to safely undergo MRI examination or allowing the patient or an individual to enter the MRI environment is 6 (six) weeks.

CLINICAL EVALUATION OF BONE PLATES & BONE SCREWS

The Miraclus Orthotech Pvt. Ltd. Bone plates & Bone Screws are clinically safe, and effective in single use as discussed and proved up to the mark in the clinical evaluation of the device.




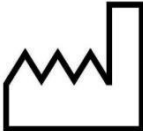



DISPOSAL OF BONE PLATES & BONE SCREWS

Please note that using a single use device (SUD) which comes into contact with human blood or tissue constitutes, these device may be a potential biohazard and should be handled in accordance with accepted medical practice and applicable local and national requirements

FOR FURTHER INFORMATION

Please contact Miraclus Orthotech Pvt. Ltd. in case of any Query, Complain or Adverse Effect

Email: info@miraclus.com / complaints@miraclus.com , Tel +91 79 4891 9551

	<p style="text-align: center;">Non-Sterile</p> <p style="text-align: center;">Indicating that the device has not been sterilized.</p>
	<p style="text-align: center;">Consult Instructions For Use</p> <p style="text-align: center;">Note: This symbol advises the reader to consult the operating instructions for information needed for the proper use of the device.</p>
	<p style="text-align: center;">Do not re-use</p> <p style="text-align: center;">Single use or use only once</p>
	<p style="text-align: center;">Date Of Manufacture</p> <p style="text-align: center;">Note: This symbol is accompanied by the date that the device was manufactured. The date could be year, year and month, or year, month and day, as appropriate.</p>
	<p style="text-align: center;">Catalogue Number</p> <p style="text-align: center;">Note: This symbol be accompanied by the catalogue number relevant to the device bearing the symbol.</p>
	<p style="text-align: center;">Batch Code</p> <p style="text-align: center;">Note: This symbol should be accompanied by the batch code relevant to the device bearing the symbol.</p>
	<p style="text-align: center;">Do Not Use If Package Is Damaged</p> <p style="text-align: center;">Do not use, if the packaging is compromised.</p>

	<p style="text-align: center;">Caution</p> <p>This symbol is to denote that there some warning or precautions associated with device, which are not otherwise found on labels</p>
<p style="text-align: center;">Qty</p>	<p style="text-align: center;">In Single Pack Number Of Quantity Packed</p>
<p style="text-align: center;">Material</p>	<p style="text-align: center;">Raw Material used for manufacturing</p>
	<p style="text-align: center;">Manufacturers Company Logo</p>
	<p style="text-align: center;">Authorized Representative in the European Community</p> <p style="text-align: center;">CMC Medical Devices & Drugs S.L</p> <p style="text-align: center;">C/Horacio Lengo N° 18, CP 29006, Malaga, Spain</p> <p style="text-align: center;">Tel: +34 951 214 054</p> <p style="text-align: center;">Email: info@cmcmedicaldevices.com</p>
	<p style="text-align: center;">Manufacturer</p> <p style="text-align: center;">Miraclus Orthotech Pvt. Ltd</p> <p style="text-align: center;">Survey No 106, Village Vavdi, Opp Safex Industries, Behind Elite Boarding School, Vavdi, Ta. Kheda-387550, Gujrat, India</p>
	<p style="text-align: center;">Medical device Indicates the item is a medical device</p>
	<p style="text-align: center;">Rx is the symbol for a medical prescription.</p>
	<p style="text-align: center;">CE marking with Notified Body Number (This symbol shall be used only after completion of conformity assessment & availability of valid CE certificate)</p>