

Basic Instrument Set

Sr. No.	Description	Code
1	Tap	01319-N
2	Drill bit 2.5 mm caliberated stopper	01320-N
3	Drill bit 2.5 mm caliberated stopper	01321-N
4	Drill bit 2.5 mm	01322-N
5	Hexagonal driver shaft	01323-N
6	Hexagonal driver shaft	01324-N
7	Cervical plate holding forcep	01325-N
8	Locking spanner for secure lock	01326-N
9	Screw holding sleeves	01327-N
10	Curette	01328-N
11	Plate bender	01329-N
12	Pin distractor	01330-N
13	Starter Awl with adjusting stopper	01331-N
14	Temporary pin - Short	01332-N
15	Temporary pin - Long	01333-N
16	Staight AO QC handle	01334-N
17	Drill guide for 2.5 mm drill bit	01335-N
18	Temporary pin holder for short pins	01336-N
19	Temporary pin holder for Long pins	01337-N
20	Plate box	-

INSTRUMENT SET



Miraclus Orthotech Pvt. Ltd.
Corporate office :
B-513, Mondeal Square, Opp. Karnavati Club,
Pralahadhnagar, Ahmedabad-380015

Plant address :
Plot No. 1112, Phase-3, G.I.D.C., Vatva, Ahmedabad - 382 445.
Phone : +91 79 2583 1332, 2583 1333
Email : contact@miraclus.com | www.miraclus.com



ACP™ System
Anterior Cervical Plating System.

- Integrated Locking Mechanism
- Increased Angulation
- Simplified Instrumentation

• Miraclus ACP System is designed to offer simplicity, speed & desired patient outcome to the spine surgeons.

• ACP System Plate is one of the slimmest profile on the market. It features an unique locking mechanism which minimises the risk of screw back out.



Adequate window for graft visibility



Integrated Locking Mechanism



Integrated locking mechanism enable with ¼ twist of a driver prevents screw backing out.

1.9 mm Uniform Thickness




Pre-contoured plate to fit the anatomy of the patient which reduces the surgery time.

Additional Bending can be possible.

18 mm

ACP Plates are available from one-level, two-level & three-level constructs.

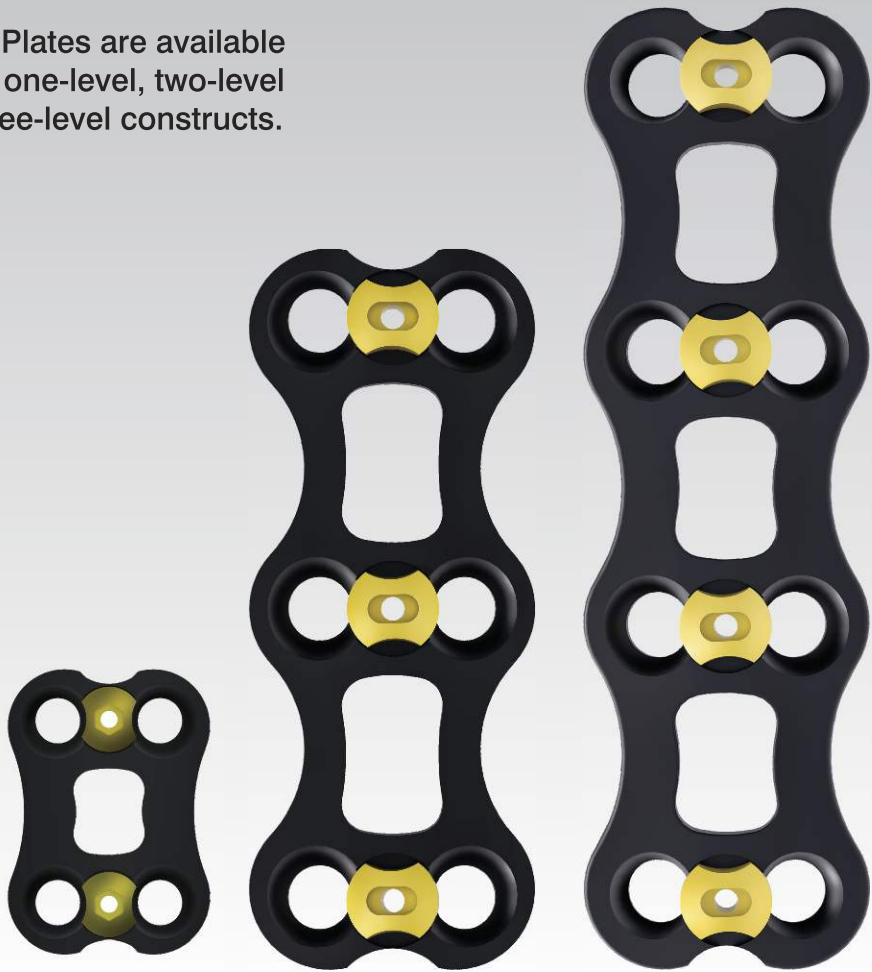

Sizes :

Material	Level	Height (mm)
Titanium Alloy	I	19, 21, 23, 25, 27.5, 30
	II	32.5, 35, 37.5, 40, 42.5, 45, 47.5, 50
	III	52.5, 55, 57.5, 60, 62.5, 65, 67.5, 70, 72.5

4.0 mm



Green

4.5 mm



Colour Coded Screw

- 4.0 mm Screws are Self-tapping, Self-drilling & green in colour for easy identification
- 4.5 mm Revision Screws are Self-tapping, Self drilling and purple in colour for easy identification.

Self- Drilling, Self –Tapping Screws

Self-Drilling, Self-Tapping Screws reduces the need for intra-operative surgical steps.

Sizes :

Description	Height (mm)
4.0 mm Cervical Screw	11, 13, 15, 17
4.5 mm Cervical Revision Screw	11, 13, 15, 17